

Innovative accessories for simplified treatment procedures

OptraGate® OptraDam®Plus OptraSculpt®Pad OptraSculpt® OptraPol® OptraFine® OptraStick®



Simplified isolation

relative isolation

OptraGate[®]

Simplified access to the cavity

- Even retraction of the lips and cheeks
- Enlarged and easily accessible treatment field
- Effective and fast relative isolation

Stress-free treatment

- Easy placement
- More time to focus on the actual treatment procedure
- Relieves the practice team and dentist

High patient comfort

- Soft and flexible material with recesses to accommodate the lip frenula
- Gentle to dental tissue also suitable for children
- Latex-free well tolerated



Source: Ivoclar Vivadent AG

Range of indications





Restorative therapy



Intraoral scans

Additional indications

ndications

- Diagnostics
- Tooth whitening
- Orthodontics
- Pediatric dentistryProsthetics



User video available at www.ivoclarvivadent.com



NEW: now available in blue and pink

absolute isolation

OptraDam[®] Plus

Good overview

- Optimum adaptation to the anatomy of the oral cavity thanks to the three-dimensional shape
- Even retraction of the lips and cheeks in all directions
- Enlarged treatment field

Easy and fast placement

- Integrated frame and perforation pattern
- Allows hold without clamps
- Contact points are easily overcome

High patient comfort

- Provides gentle retraction and protection of the soft tissues
- Low coverage of extraoral tissues









Simplified contouring

anterior region

OptraSculpt[®]**Pad**

Effortless adaptation due to anti-stick effect

- Reduced adhesion of the composite to the instrument
- Efficient and easy contouring

Professional esthetics

- No instrument marks on the restoration surface
- Allows for a homogeneous surface design
- Additional reference scales on the instrument handle provide orientation for the axial position and width of anterior teeth
- Suitable for the build-up of direct veneers or large Class IV restorations

High efficiency

- Faster contouring
- Reduced effort required for finishing and polishing





Source: Dr Ronny Watzke, Liechtenstein, 2014





User video available at http://www.ivoclarvivadent.com/embed/16237

posterior region

OptraSculpt[®]

Simplified adaptation due to non-stick effect

- Reduced adhesion of the composite to the instrument
- Efficient and convenient contouring

High flexibility

- 6 attachments for different indications
- Modelling tips can be rotated by 360° and snapped into place at any angle desired

Efficiency increase

- Faster contouring
- Reduced effort required for finishing and polishing



Simplified polishing

composite

OptraPol[®]

1-step polishing system


- Radiant high-gloss restoration surfaces in only one step
- High diamond content (approx. 72 wt.%) for excellent results

High efficiency

- High gloss after a short time and in only one step
- No need for changing tips

Long-lasting and economical

- High dimensional stability
- High wear resistance
- Autoclavable and reusable
- Cost-saving
- 4 shapes for different treatment areas

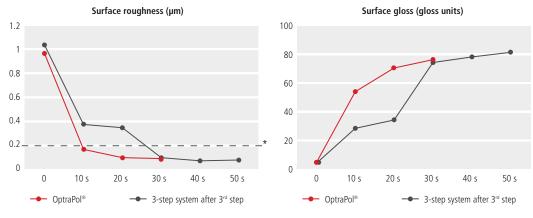

* clinically relevant threshold value (Bollen, C. M. et al, 1997)

Sources:

In vitro test of competitors' polishing systems; composite: Tetric EvoCeram® (Dr S. Heintze, Ivoclar Vivadent, Schaan, 2009)

Bollen, C. M., Lambrechts, P. & Quirynen, M. (1997). Comparison of surface roughness of oral hard materials to the threshold surface roughness for bacterial plaque retention: a review of the literature. Dent Mater, 13(4), 258-269.



ceramics

OptraFine[®]

Highly coordinated with ceramic materials

- High-performance polishing system
- Suitable for ceramics (e.g. including IPS Empress[®] CAD and IPS e.max[®] CAD)
- Particularly smooth restoration surfaces

3-step polishing system

- Step 1: Pre-polishing using "light blue"
- Step 2: Polishing using "dark blue"
- Step 3: High gloss with brushes and polishing paste

Long-lasting and economical

- High wear resistance
- Autoclavable and reusable
- Cost-saving
- 3 shapes for different treatment areas



Simplified placement

indirect restorations

OptraStick[®]

Precise positioning

Application

- Adhesive tip allows accurate placement of restorations
- Easy pick-up and release
- Particularly suitable for the positioning of small indirect restorations

Suitable for almost any material

- Inlays, onlays, crowns and veneers made of ceramic (e.g. IPS Empress[®] or IPS e.max[®])
- Veneers made of IPS d.SIGN[®] and IPS InLine[®]
- Inlays, onlays and ³/₄ crowns made of gold



Delivery forms



OptraGate Assortment 40 x size R (regular) 40 x size S (small)

Refill R 80 x size R (regular) Refill S 80 x size S (small) Refill Junior 80 x size Junior



OptraGate in blue and pink

Junior Assortment/blue and pink 20 x size Junior blue

20 x size Junior pink

Starter Kit/blue and pink

4 x size Junior blue 4 x size Junior pink 4 x size S (small) blue 4 x size S (small) pink

Small Assortment/blue and pink 20 x size S (small) blue

20 x size S (small) pink



OptraDam Plus

Assortment 25 x size R (regular) 25 x size S (small)

Refill R

50 x size R (regular) **Refill S** 50 x size S (small)



OptraSculpt Pad Assortment

1 x OptraSculpt Pad Instrument 1 x OptraSculpt Pad Refill/4 mm/60 1 x OptraSculpt Pad Refill/6 mm/60

OptraSculpt Pad Refill/4 mm/60 OptraSculpt Pad Refill/6 mm/60 OptraSculpt Pad Instrument Refill



OptraSculpt Assortment

1 x Modelling instrument (single-ended) 1 x Modelling instrument (double-ended) 40 x Modelling tips (ball) 40 x Modelling tips (point) 40 x Modelling tips (spatula) 40 x Modelling tips (chisel) 40 x Modelling tips (cylinder) 40 x Modelling tips (pyramid)

Refill Modelling Instrument

1 x Modelling instrument (single-ended) **Refill Modelling Instrument** 1 x Modelling instrument (double-ended) **Refill Modelling Tips** Each shape is separately available in packets of 100.

OptraPol

Assortment

- 6 x Polishers small flame
- 4 x Polishers large flame
- 4 x Polishers cup
- 4 x Polishers lens

Refill Polisher

Each shape is separately available in packets of 10.

OptraFine

Assortment

- 2 x Flame F (finisher) 2 x Flame P (polisher)
- 2 x Flame I
- 2 x Cup F
- 2 x Cup P
- 2 x Disc F
- 2 x Disc P
- 6 x Nylon brushes
- 1 x Polishing paste

Refills

Each polishing instrument is separately available in packets of 10. Polishing paste refills are available in single syringes.

OptraStick

Refills of 50



Direct Restoratives

Optra Accessories form a part of the "Direct Restoratives" product category. The products of this category cover the procedure involved in the direct restoration of teeth – from preparation to restoration care. The products are optimally coordinated with each other and enable successful processing and application.



THESE ARE FURTHER PRODUCTS OF THIS CATEGORY:

Tetric EvoCeram[®] Bulk Fill & Tetric EvoFlow[®] Bulk Fill

High-performance posterior composite



The most efficient posterior composite!

- Increment thickness of up to 4 mm due to the highly reactive light initiator lvocerin $^{\circ}$
- Dentin- or enamel-like volume replacement
- 10 seconds (≥ 1,000 mW/cm²)
- 47% less time required than with the conventional technique* * compared to Tetric EvoFlow® and Tetric EvoCeram®. Data available on request.

Bluephase® Style

The curing light



The smallest LED for every use

- Every hand comfortable to hold for men and women
- Every material universal use due to polywave® LED with broadband spectrum
- Every time optional corded operation due to Click & Cure

Would you like to know more about the products of the "Direct Restoratives" category? Simply get in touch with your contact person at Ivoclar Vivadent or visit www.ivoclarvivadent.com for more information.

Ivoclar Vivadent AG Bendererstr. 2 9494 Schaan Liechtenstein Tel. +423 235 35 35 Fax +423 235 33 60 www.ivoclarvivadent.com

Descriptions and data constitute no warranty of attributes. Printed in Liechtenstein © Ivodar Vkadent AG, Schaan, Liechtenstein. 640990/2015-06/8/KNRev, 1

